
TERMOMETRU DIGITAL DIGITAL THERMOMETER

TERMOMETRU DIGITAL

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DIGITAL THERMOMETER

Thank you for purchasing the Trussa Care for all digital thermometer. This is a medical device for measuring human body temperature.

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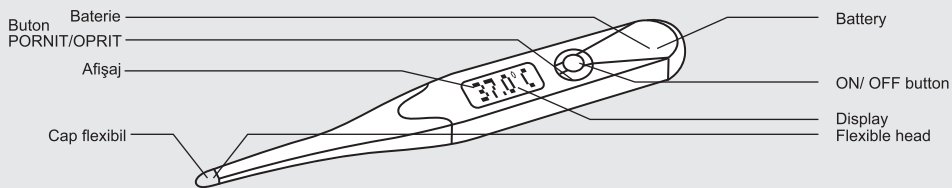Directions for use:

1. Press the ON/OFF button to start the thermometer.
2. Place the head of the thermometer in the desired area. Wait for the sound signal.
3. Read the result.
4. Press the ON/OFF button to stop the thermometer.

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METHODS TO MEASURE THE TEMPERATURE

ORAL USE

Clean the flexible head.
Place the flexible head under the tongue.
Keep the flexible head under the tongue for about 30 - 60 seconds, until you hear the sound signal.
The result on the display is final.

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UNDERARM (AXILLARY) USE

Clean the flexible head.
Place the flexible head in the armpit.
Keep the flexible head in the armpit for about 1-2 minutes, until you hear the sound signal.
Add 0.5°C to the result displayed, in order to read the final result.

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RECTAL USE

Clean the flexible head.
Place the flexible head into the rectum.
Keep it in the rectum for about 1 minute, until you hear the sound signal.
Subtract 0.5°C from the result displayed, in order to read the final result.

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CLEANING THE THERMOMETER

You can disinfect the head of the thermometer using water and soap or rubbing alcohol any time you use it for measuring the temperature.

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ATTENTION!

1. Keep away from chemical solvents, direct sunlight or high temperatures.
2. Do not recharge the battery and do not throw it in fire.
3. For hygienic purposes, clean the head of the thermometer before every use.
4. The temperature should be measured only under adult supervision.

PRODUCĂTOR/ MANUFACTURER: ESTAR ELECTRONIC CO., LTD., Xindu Industry Zone, Changshan Country 324200, China
pentru for: PHARMATECH ASIA GROUP LTD., 8A, City Hotel Building, Minle Industry Park, Meiban Road, Longhua New District, Shenzhen, China 518131.

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